

# GEM MicroClip

**GEM Precision Tools**  
for the enhancement of microsurgery

## Have you ever had a hemostatic clip fall off?

The GEM titanium, hemostatic MicroClip (GEM2431) is 30% smaller than the small-sized clips presently on the market. The secure closure of the MicroClip is formed by its Chevron shape. The tip-to-tip closure is progressive and encompassing, and the atraumatic triangular design of the MicroClip ensures a uniform closure. The inside of the MicroClips are lined with diamond-shaped grooving that create interlocking teeth that enhance security and vessel closure while preventing multi-directional slippage. GEM MicroClips are great for side branch occlusions.

## Benefits of the GEM MicroClip:


- The MicroClips are stored in a precise and simple cartridge that ensures easy access.
- The MicroClips have an optimal, atraumatic clip design.
- The MicroClips are pure, strong, malleable titanium.
- The MicroClips are MRI compatible.

## The GEM MicroClip System consists of:

GEM2431      MicroClips; 30 cartridges of 6 per box

GEM615-1G    MicroClip Applier; 15 cm length

GEM619-1G    MicroClip Applier; 19 cm length

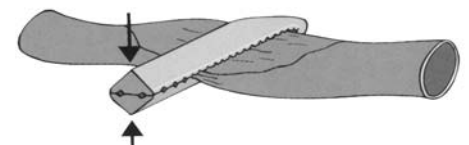

Ultra-smooth, atraumatic, triangular-formed wire ensures uniform clip closure



Chevron shape provides progressive and encompassing tip-to-tip closure before complete and secure occlusion of tissue or vessel



Diamond-shaped grooving creates atraumatic, interlocking teeth that enhance security and vessel closure while preventing multi-directional slippage

The GEM MicroClip is manufactured by Vitalitec International and distributed by Synovis Micro Companies Alliance, Inc.

To place an order or for technical support, please contact Synovis MCA at 1-800-510-3318.

Visit our website at [www.synovismicro.com](http://www.synovismicro.com)

**Synovis®**

Micro Companies Alliance, Inc.

439 Industrial Lane

Birmingham, AL 35211-4464 USA

Tel: 205.941.0111 or 1.800.510.3318

Fax: 205.941.1522

[synovismicro.com](http://synovismicro.com)