

Orange Means Micro MICROTRAY™

Protecting Your Investment The GEM MicroTray Organizes & Protects Microsurgical Instruments During Sterilization & Storage

Anodized Aluminum Base And Cover

- Light weight aluminum for strength and impact resistance
- Anodized and coated to provide an attractive, durable finish
- Silicone inserts and retaining strips are Latex-Free

Perforated Tray Base, Lid, And Inserts

- Facilitates steam penetration and moisture removal during sterilization
- Intended to be used with a sterilization wrapper to maintain sterility, or
- Fits easily into most rigid container systems
 - Single layer tray dimensions: L20.5" x W10.25" x H2"
 - Double layer tray dimensions: L20.5" x W10.25" x H3.5"

Protecting Your Investment

- Silicone finger mats and retaining strip on the inside of the lid ensures secure placement during transport
- Available with silicone instrument rack for the GEM Microvascular Anastomotic COUPLER System and GEM MicroClip and SuperFine Clip Appliers

Backed by a Two Year Warranty

- All metal parts are warranted for two years from the date of purchase against any defects in function or usability
- Refer to the Instructions For Use or Package Insert for Care and Handling Instructions

Ordering Information

GEM2551	Single layer tray	Full finger mat
GEM2551-C	Single layer tray	Finger mat and instrument rack
GEM2551-2	Double layer tray	Full finger mats upper and lower tier
GEM2551-2C	Double layer tray	Upper tier finger mat, lower tier finger mat and instrument rack



Synovis®


Micro Companies Alliance, Inc.

Instrument Repair Program

In need of repair? Trust your little GEMs with Synovis Micro Companies Alliance, Inc. This will ensure that repairs are performed properly using only factory authorized parts and certified personnel.

1. Call **800-510-3318** to obtain an RGA #
2. The instrument must be clean, sterile and in a sterile peel pouch with indicator. All items will be sent back if proper sterility is not visible.
3. The **COUPLER** instrument must be manufactured by Synovis Micro Companies Alliance, Inc.
4. **S&T** instruments must be sent in separate sterile peel pouches with indicator. Do not send instruments together as they may become damaged during transport.
5. The instrument needs to be sent alone – the tray and other components should not be sent.
6. After receiving and entering in the RGA#, ship to:
Attn: Repairs /RGA # _____
Synovis Micro Companies Alliance, Inc.
439 Industrial Lane
Birmingham, AL 35211
7. Please reference your RGA # on the outside of the shipping container.

*Repairs take approximately 2 – 3 weeks.
If you need a loaner, let our customer service representative know.*

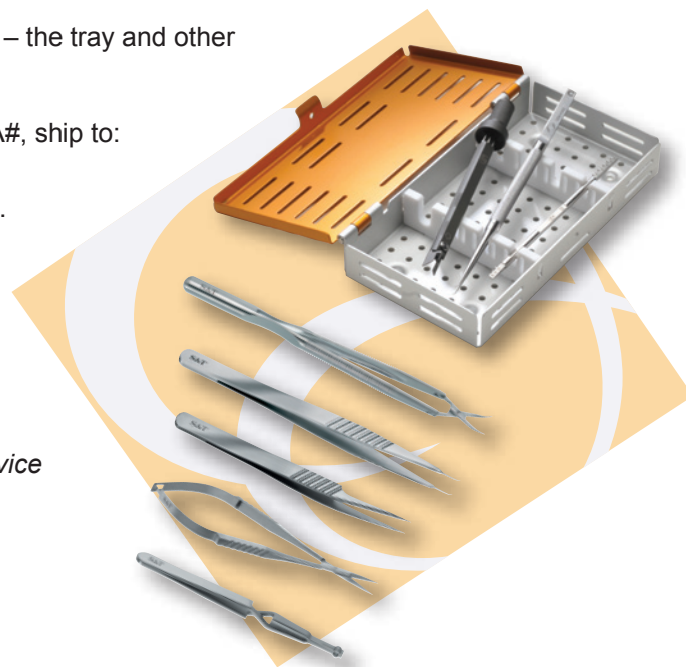

Manufactured for: Synovis Micro Companies Alliance, Inc.,
a Subsidiary of Synovis Life Technologies

